

Pipeline for Comparison of Weighted Structural Data between Groups of Subjects

The BRAPH 2 Team

September 29, 2023

In this tutorial, we will upload a file containing the pipeline with the different steps to compare two groups of subjects using *structural data* (check tutorial *Group of Subjects with Structural Data*) and weighted undirected graphs. This Tutorial explains how to perform group analyses and comparisons with this kind of data.

Contents

<i>Open the GUI</i>	3
<i>Step 1: Load the Brain Atlas</i>	3
<i>Step 2: Load the Structural Group Data</i>	4
<i>Step 3: Analysing the Data of Group 1</i>	4
<i>Step 4: Analysing the Data of Group 2</i>	7
<i>Step 5: Comparing Groups</i>	7

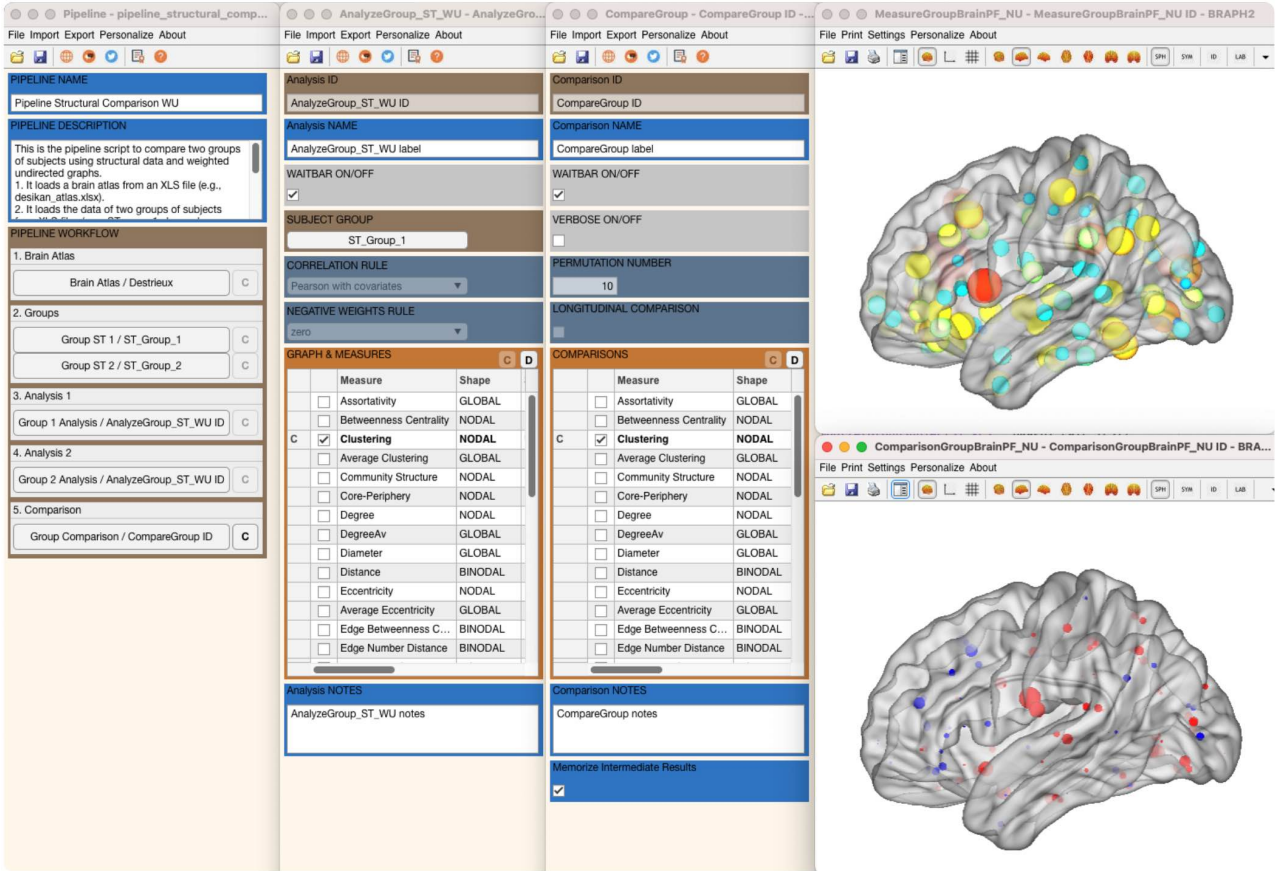


Figure 1: GUI for working with the pipeline to compare two groups of subjects with *structural data* using a weighted undirected analysis (Figure TBD). Full graphical user interface to perform group comparisons with *structural data* in BRAPH 2.0.

Open the GUI

The pipeline GUI can be opened by typing `braph2` in MatLab's terminal, which allows you to select a pipeline, in this case, `pipeline_structural_comparison_wu.braph2` (containing the steps required to perform your analysis).

To open the GUI and upload the structural comparison pipeline, you can also do it from the command line by typing the commands in Code 1.

Code 1: Code to launch the GUI to upload a pipeline file to compare two groups of subjects. This code can be used in the MatLab command line to launch the GUI to upload a pipeline file.

```
1 pip = Pipeline(); ①
2
3 gui = GUIElement('PE', pip); ②
4 gui.get('DRAW') ③
5 gui.get('SHOW') ④
```

① creates a new object `Pipeline` to upload the steps to perform an analysis or comparison.

② creates a GUI to upload the files into the pipeline.

③ draws the GUI.

④ shows the GUI.

Once the pipeline is uploaded, you can see a GUI that contains different steps: to upload a brain atlas, to upload the structural data of two groups, analyze them, and finally, to compare the groups (Figure 2a).

Step 1: Load the Brain Atlas

Figure 2 shows how to upload and plot the brain atlas that you used to extract the *structural data* for your analysis. For more information on where to find different atlases or how to change plotting settings on the brain surface, check the *Brain Atlas* tutorial.

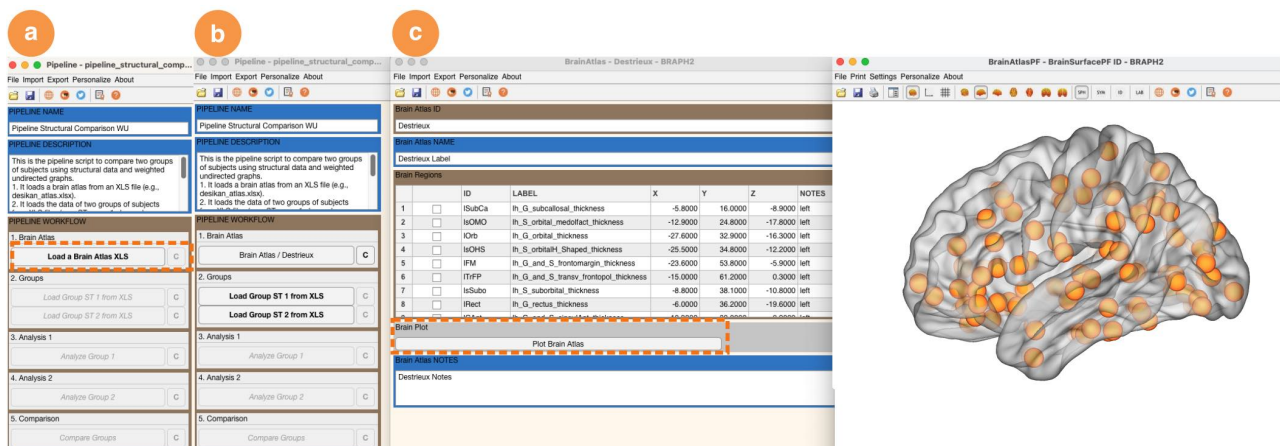


Figure 2: Uploading the Brain Atlas. Steps to upload the brain atlas: **a** Click on Load Atlas from the pipeline GUI. **b** Navigate to the BRAPH 2.0 folder atlases and select one of the atlas files, in this example the `destrieux_atlas.xlsx`. **c** You can also plot the brain atlas by pressing Plot Brain Atlas.

Step 2: Load the Structural Group Data

After you loaded the brain atlas, you can upload the *structural data* for each group as shown in Figure 3. A new interface will be shown containing the data for the group you just selected. You can open each subject's structural values by selecting the subject, right click, and select "Open selection" (for more information check tutorial *Group of Subjects with Structural Data*).

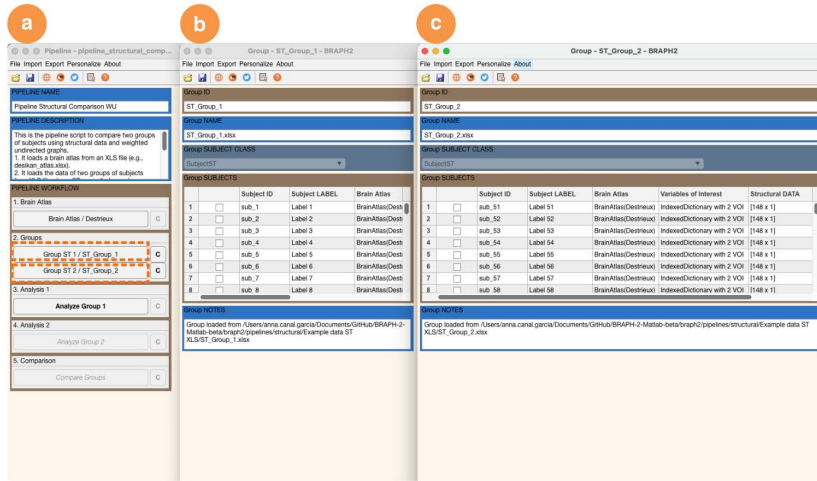


Figure 3: **Loading the Groups Data.** **a** From the pipeline GUI, click on Load Group ST 1 XLS to load the data of group 1, and Load Group ST 2 XLS to load the data of group 2. **b** Data for group 1 is uploaded. **c** Data for group 2 is uploaded.

Step 3: Analysing the Data of Group 1

After you have loaded the data for both groups, you can analyze the data of your first group by clicking on Analyze Group 1 (Figure 4a-b). A new interface will be shown that allows you to pre-calculate network measures for each group and explore them. First of all, you can specify the parameters for constructing the graph from the structural group data (Figure 4c): you can select the statistical test to calculate correlations in structural values between pairs of brain regions (Figure 4d), and you can decide how you want to analyze the negative weights from the correlation (Figure 4e).

Importantly, the parameters you select for this analysis (graph and measures parameters) will be fixed for the analysis of the second group and for the group comparisons.

After setting these parameters, you can calculate some graph measures (Figure 5). First of all, you need to ensure the measure is using the rule or parameter that you wish (not all measures have rules or parameters). As an example, we can select the measure Clustering and right-click on the top of the table (Figure 5a), in the GRAPH & MEASURES

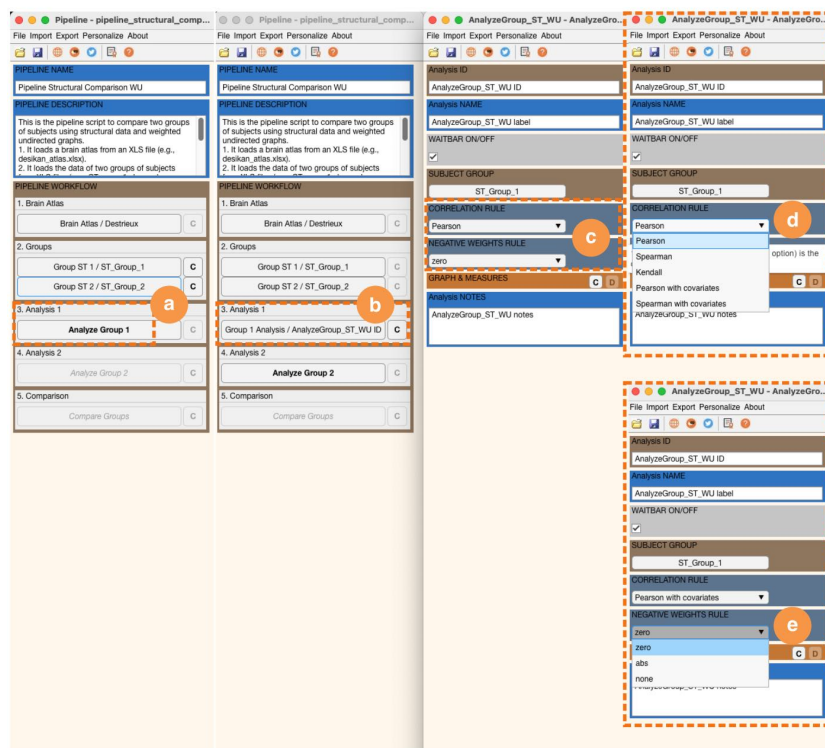


Figure 4: **Analyze the Group Data.** **a** Click on Analyze Group 1 in the pipeline's GUI. **b** After clicking on Analyze Group 1 in the pipeline's GUI, a new window will appear where first you can select the parameters for the graph construction **c**. In this pipeline, you can select the statistical test to be used for the correlations **d**, and what you want to do with the negative correlation weights **e**.

panel, and press “Data Selected Measures”. Next, the measure window will appear, allowing you to choose the rules or parameters (Figure 5b). In this instance, the rule applies to the calculation of triangles in directed graphs. However, since we’re working with an undirected graph, we will keep the default option (Figure 5b). After the rule is set, you can calculate the measure by pressing the “C” (Figure 5c) and the results will appear in a new table within the same panel (Figure 5d). Also notice that after the measure is calculated, the rule is blocked (Figure 5d).

If you want to visualize the results, select the measure and press Plot Selected Measures on Brain in the analysis’ GUI (Figure 6a). In settings (Figure 6b-c), you can change the visualization of the plots and save them.

Finally, when you do right click, in the GRAPH & MEASURES panel, there are other options you can explore such as Plot Graph Plot (connectivity adjacency matrix) as well as Data Graph (labels of brain regions, values of the adjacency matrix, options to plot matrix and histogram), all of which can also be saved.

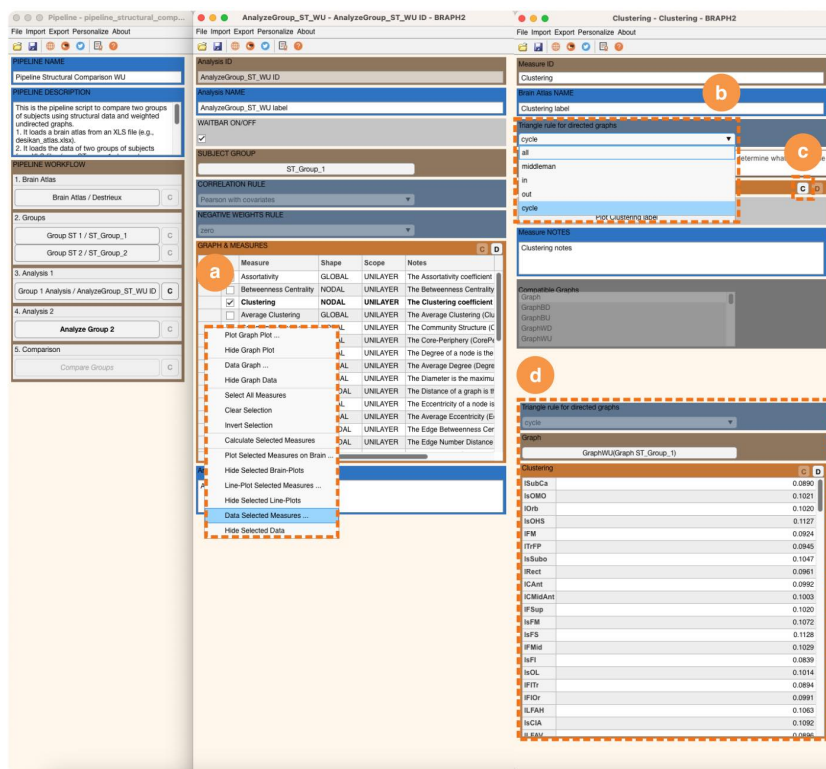


Figure 5: Check and set the rules and/or parameters of a measure before its calculation. **a** Click on Data Selected Measures after selecting the measure. **b** This new window shows the array of values for the measure, in this case Clustering. **c** This new window shows the array of values for the measure, in this case Clustering. **d** After the calculation of the measure, a new table appears with the measure results, and the rule is blocked (in blue).

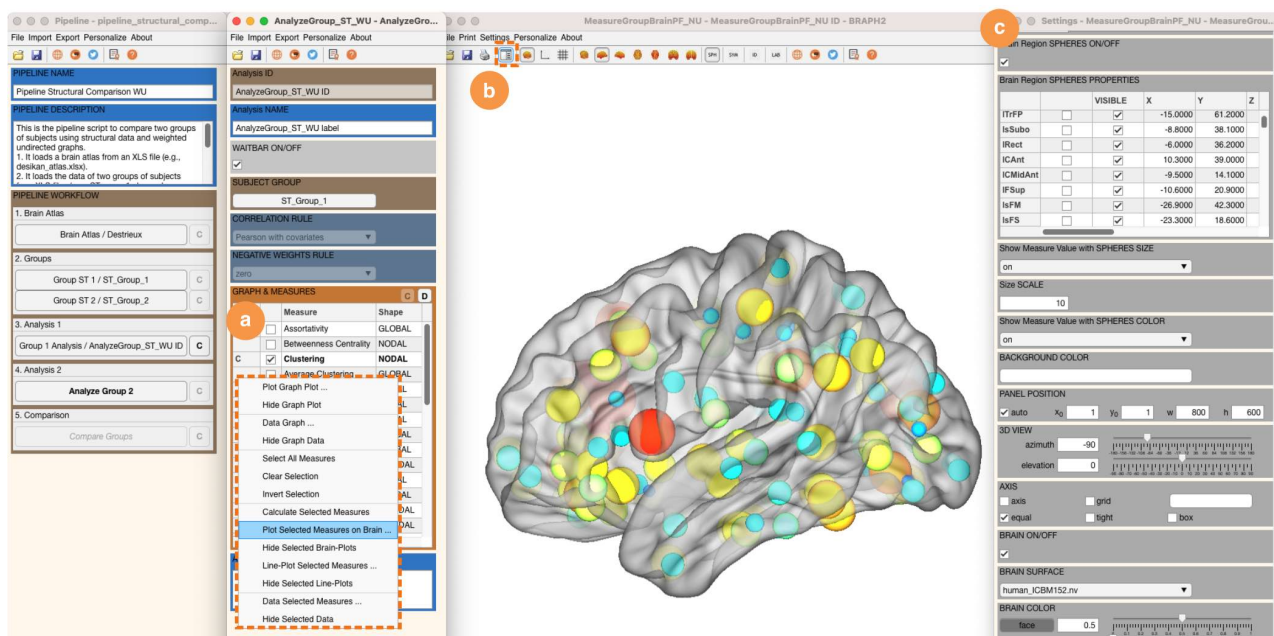


Figure 6: Visualize the measure's results in a brain plot. **a** Click on Plot Selected Measures on Brain in the analysis' GUI and a brain surface with the results from the calculated measures will appear. The size of the spheres and the color are proportional to the measure's value. **b** In the new window, press the button settings to obtain further visualization options of the results **c**.

Step 4: Analysing the Data of Group 2

After the analysis of group 1, you can proceed with the analysis of the second group by clicking on Analyze Group 2 (Figure 7a). You will notice that, in the new window that is shown (Figure 7b-d), the parameters you selected for group 1 are already selected and fixed for this analysis (both graph and measure parameters). If you realize that some of the options you previously selected are not the ones you would like, you can reset the analysis parameters of group 1 by clicking on the C checkbox next to it.

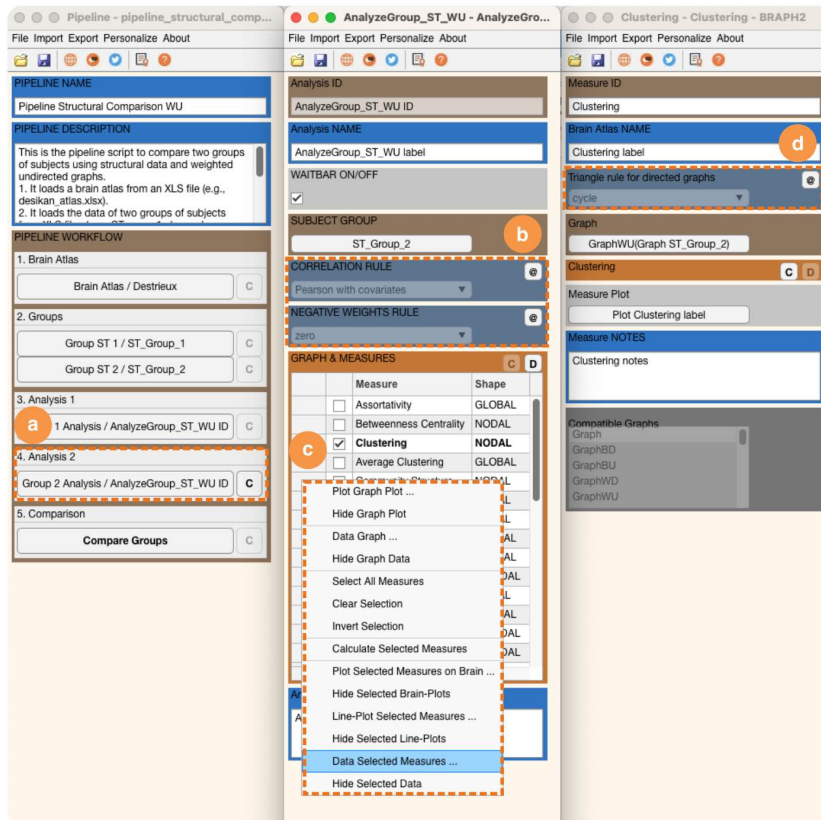


Figure 7: **Parameters blocked in Analysis of Group 2.** **a** Click on Analysis 2 in the pipeline's GUI. **b** In this new window, you can see that the graph properties such as Correlation rule and Negative weights rule are blocked since they are the same as the ones set in the analysis of group 1. If you select a measure, in this case Clustering, and press Data Selected Measures **c**, you can observe that the measure's rules and parameters are also set in case you calculated the measures in analysis 1, and if not, you can set the rule by pressing at the "@" **d**.

Step 5: Comparing Groups

Once you have explored the network measures for each group, you can proceed with their statistical comparison. To do this, you should click on Compare Groups (Figure 8a) and in the new window select if you want a waiting bar and verbose functions ON while you wait for the analysis to finish, as well as how many permutations you want to use to assess differences between groups (Figure 8b), and if the groups are not independent but represent the same subjects in two differ-

ent points in time, you can select the longitudinal comparison option, which will permute the values within each subject (Figure 8b). We set the permutations to 10 for computational time purposes (Figure 8c), but for your research analysis we recommend using 1000 or 10000 permutations. Finally, you can select the graph measures you want to compare between groups and once you have selected all the measures you are interested in, you should right click and select Calculate all selected comparisons (Figure 8d). If you turn ON the wait bar and verbose functions, two window bars will open that show you at which point in time the comparison calculation is. There is one last option on this GUI that you can select to save intermediate results during the permutations.


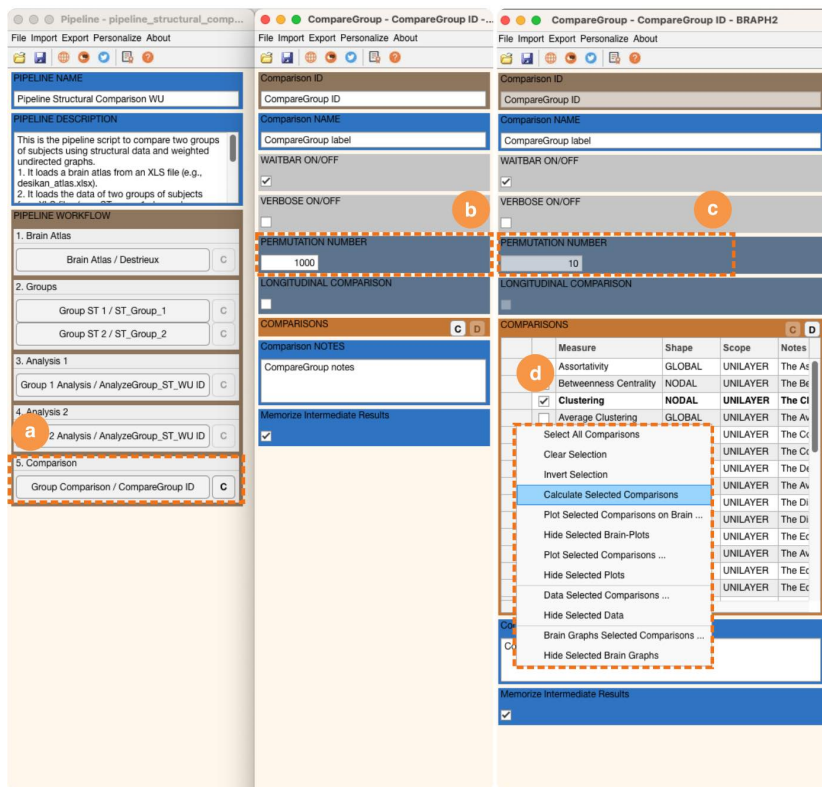


Figure 8: **Compare the Group Data.** **a** Click on Compare Groups in the pipeline's GUI. **b** In this new window, you can select what to turn ON/OFF the wait bar and verbose functions, you can change the number of permutations, and whether to perform a longitudinal group comparison. We set the number of permutations to 10 for this tutorial **c**. Finally, you can calculate the comparison of some graph measures between groups **d**.

To obtain the results from the measure/s comparison, select the measures in the GRAPH & MEASURES panel and press Data Selected Comparisons (Figure 9a), and a new window will open (Figure 9b) where we can check the difference value between groups, the p-values (1-tailed and 2-tailed), as well as the confidence intervals.

Finally, we can visualize the comparison results on a brain surface by selecting the measure comparisons we want and right click

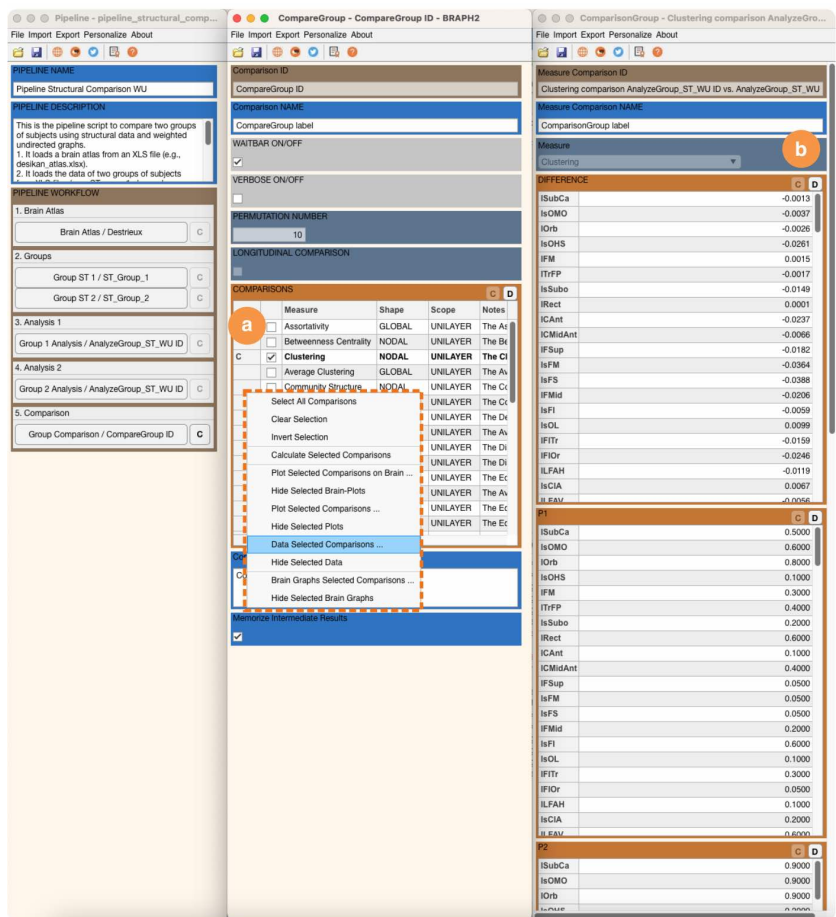

and press Plot Selected Comparisons on Brain (Figure 10a). A new window with comparison results on a brain surface will appear, where blue color indicates group 1 > group 2 and red color group 2 > group 1, and the size of the spheres is proportional to the absolute difference value between groups. If you press the Settings button (Figure 10b), then you will have available more options, for example to apply FDR correction to your plot (Figure 10c), which by default is set to off.

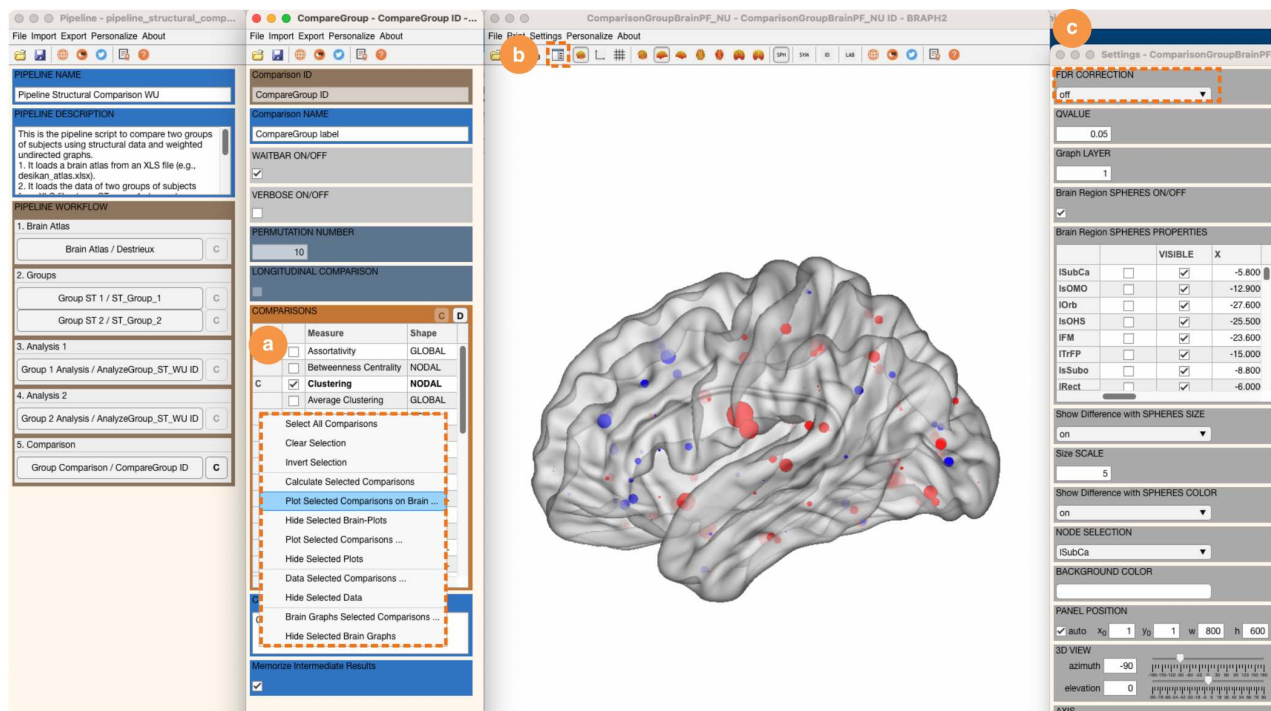


Figure 10: **Visualize the comparison results on a brain surface.** **a** Click on Plot Selected Comparisons on Brain in the comparison's GUI and a brain surface with the results from the calculated measures comparison will appear. The size of the spheres is proportional to the measure value difference between groups and the colors indicate if it is a positive or negative difference. **b** In the new window, press the button settings to obtain further visualization options of the results. **c** In the settings, you can activate the FDR correction by setting the desired q value and changing from off to on the button.